



SF40A/BT

Bluetooth Impact Printer



OHAUS Wireless Bluetooth Impact Printer

The new SF40A/BT is a portable printer that offers high performance, advanced functionality, and wireless Bluetooth connectivity to OHAUS balances and scales – making it suitable for environments that require cable-free printing, such as biological laboratory cabinets in the cleanroom.

Unique Features:

- Secure and simple setup, with a factory paired Bluetooth adapter leading to plug-and-play integration.
- With its statistics and totalization modes and clock, the SF40A/BT ensures accurate and compliant reporting. These features complement OHAUS weighing instruments' application modes and provide increased functionality for laboratory, industrial, and educational applications.
- With the ability to print in 13 different languages, SF40A/BT can print in the language of the connected device to support universal operation.



Do more with OHAUS

Check out the OHAUS laboratory equipment portfolio!

www.ohaus.com

SF40A/BT Bluetooth Impact Printer

Model	SF40A/BT
Item number	30944359
Number of Columns	40 columns
Print Speed	Approximately 1,2 lines per second
Application Modes	Totalization, Statistics (Statistics printout includes number of samples, average value, standard deviation, minimum, maximum, min-max differential, and total sample weight sum)
Paper Roll (W × Dia.)	Standard paper 57,5 mm × Ø 60 mm max.
Power Requirements	AC adapter Input: 100 – 240 VAC, 50 – 60 Hz AC adapter output: 12 VDC, 2,5 A
Interface	1 × RS232, 2 × Bluetooth Adaptors
Ambient Temperature	0 °C – 40 °C
Printer Dimensions (W × L × H)	140 × 205 × 75 mm
Shipping Dimensions (W × L × H)	255 mm × 205 mm × 150 mm
Wireless standard*	2,1
Active Communication Range*	10 m (without obstacles)
Net Weight	824 g (paper roll included)
Shipping Weight	0,8 kg

Other Standard Features and Equipment

Status indicator, 40 characters per line, perform zero and tare operations directly from the printer. Full functionality with Explorer™, Navigator™, MB Series Moisture Analyzers and 'T' Series indicators. Statistics and Totalization not supported on Scout™ Series. Multi-Language Printing: English, Spanish, German, Italian, French, Chinese, Korean, Japanese, Russian, Polish, Portuguese, Hungarian and Czech.


Compliance

- **SF40A printer:** Product Safety: IEC/EN 61010-1, : UL61010-1; CAN/CSA C22.2 61010-1
- **Electromagnetic Compatibility:** FCC Part 15 Class A; Canada ICES-003 Class A; IEC/EN 61326-1 Class B, basic environments
- **Compliance Marks:** CSA, CE
- **Bluetooth adapter certification:** Anatel (for Brazil), CE RED (for EU), FCC (for US), IC (for Canada), iDA (for Singapore), KCC (for Korea), MoC (for Israel), NBTC (for Thailand), NCC (for Taiwan China), NOM (for Mexico), QAS SIRIM (for Malaysia), RCM (for Australia and New Zealand), SRRC (for China), Telec (for Japan), Bluetooth SIG.

Accessories

Paper Roll × 2 (included)	12120799
Paper Rotary Axis (included)	30063920
Ink Ribbon (included)	30529322
Cable Adapter (included)	30059316
Cable Mini USB	30944363

Outline Dimensions

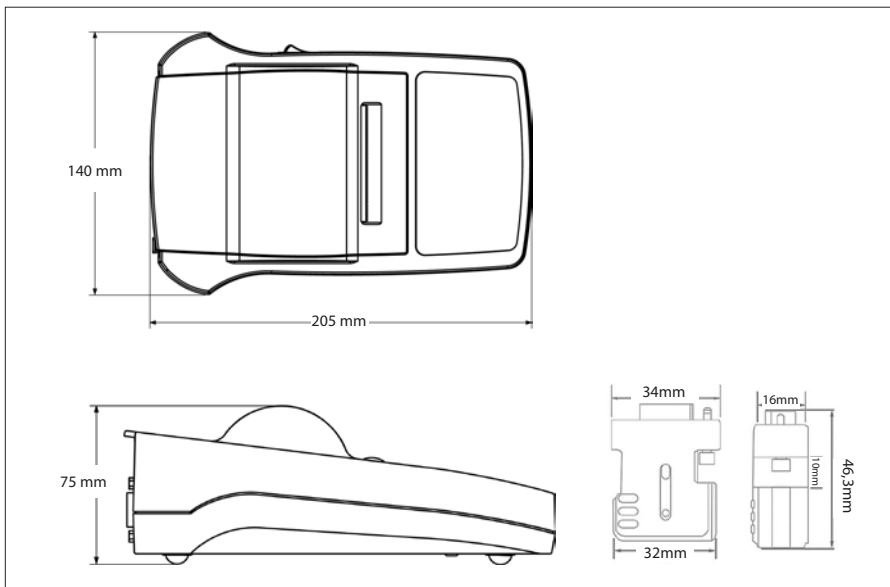

80776945_A 20240627 © Copyright OHAUS Corporation

OHAUS Europe GmbH

e-mail: ssc@ohaus.com

Tel: 0041 22 567 53 19

e-mail: tsc@ohaus.com

Tel: 0041 22 567 53 20

www.ohaus.com

The management system governing the manufacture of this product is ISO 9001:2015 certified.

